



DG5000 Series Self Test Procedure

Date:08/12/2011

Solution: To perform a self test on the DG5000 series of arbitrary waveform generators, follow these steps:

- 1. Power on the unit.
- 2. Press Utility



- 3. Right arrow





4. Press Test/Cal button




5. Press Secure Code






- 7. Enter the security code “2010”. Start by using the rotating selector wheel to scroll through the characters to select a single value.. then press the right arrow to move to the next character position. Each selected/entered value will become a “*” after pressing the right arrow.

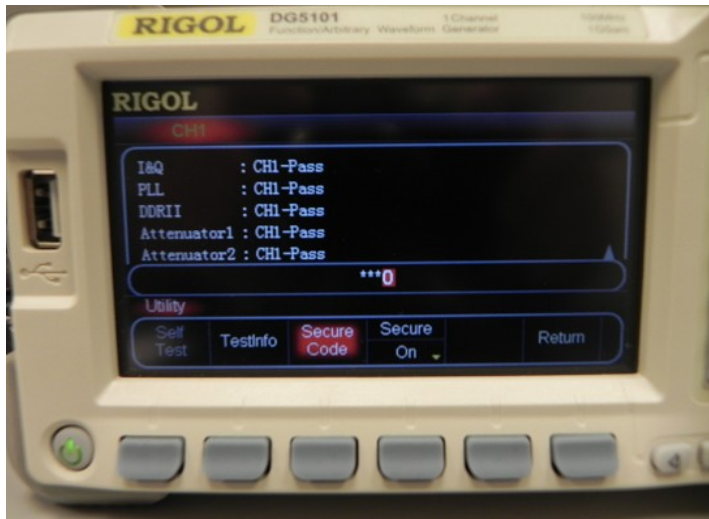



- 8. After the password has been entered, press Secure key. The unit should indicate that access has been granted.

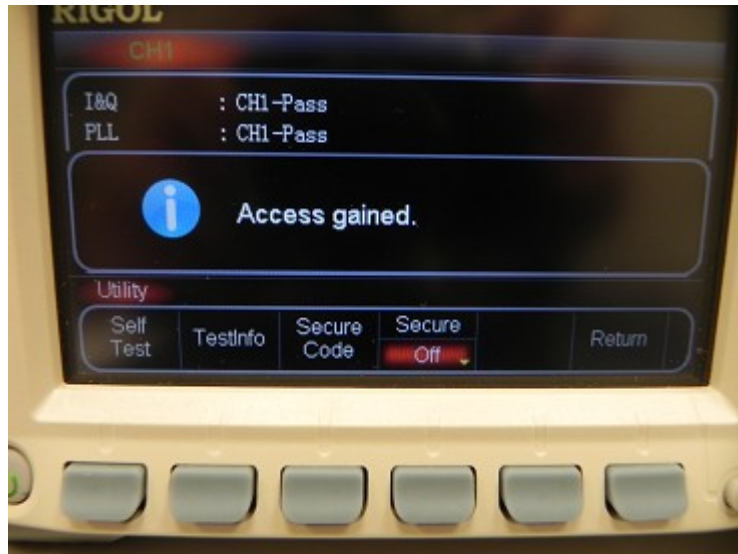

- 9. Now, you can perform a self test by pressing the Self Test button.





rigolna.com

Rigol Technologies USA, Inc.
Phone/Fax: 440-232-4488

Toll Free: 877-4-RIGOL-1

7401 First Place, Ste. N Oakwood Village, OH 44146 USA